

A full range calibration for CALIFA





Mrunmoy Jena

R3B Week
12.11.2024



Overview and structure of macro

Calibration in the gamma range



Extrapolation to proton range

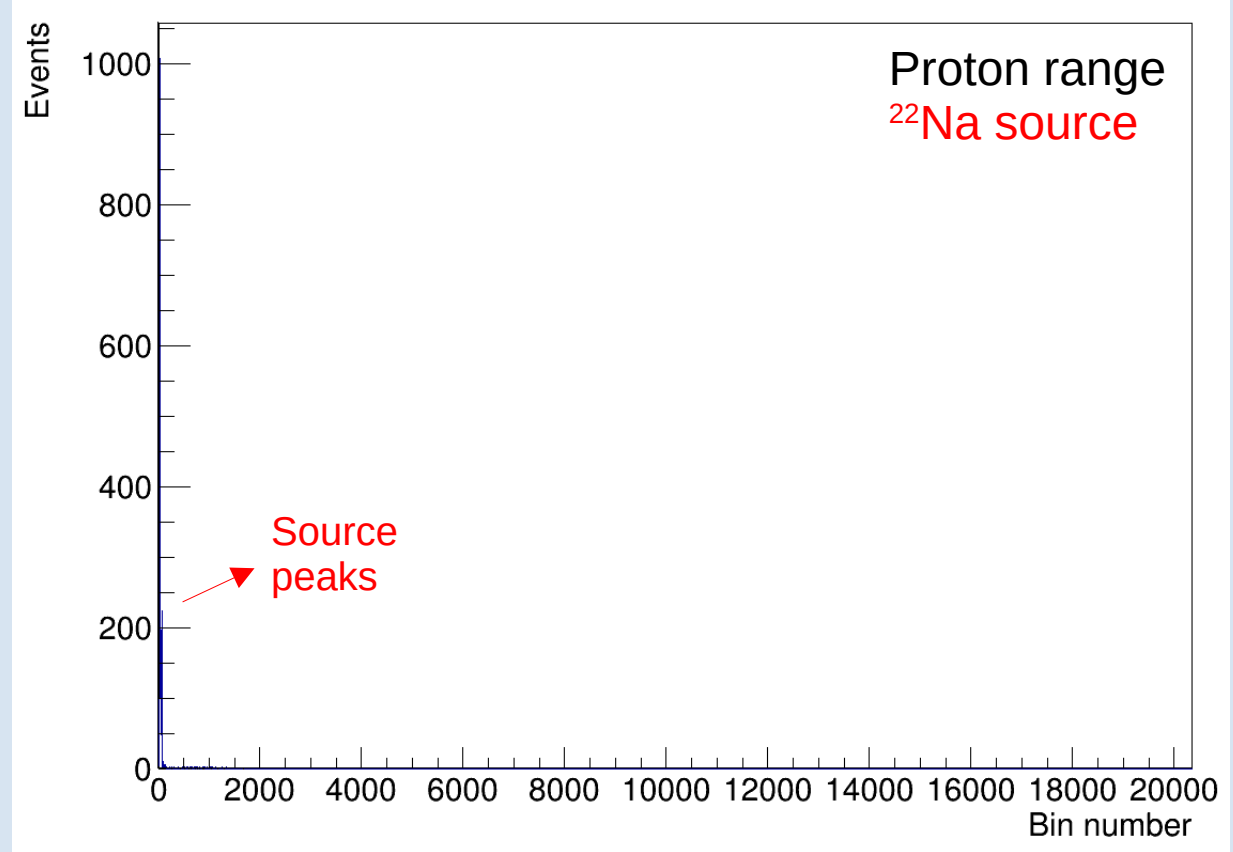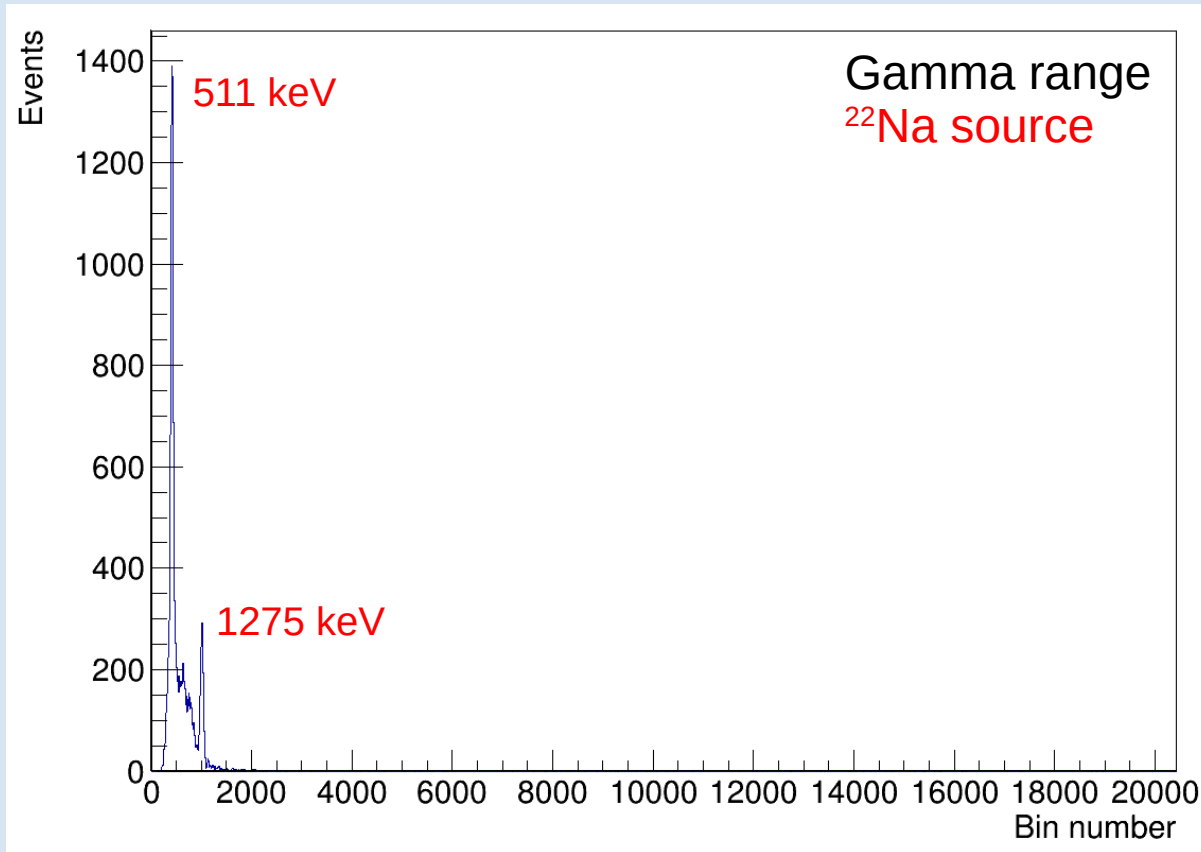


^{22}Na source hardly useful for performing calibration in the proton range !

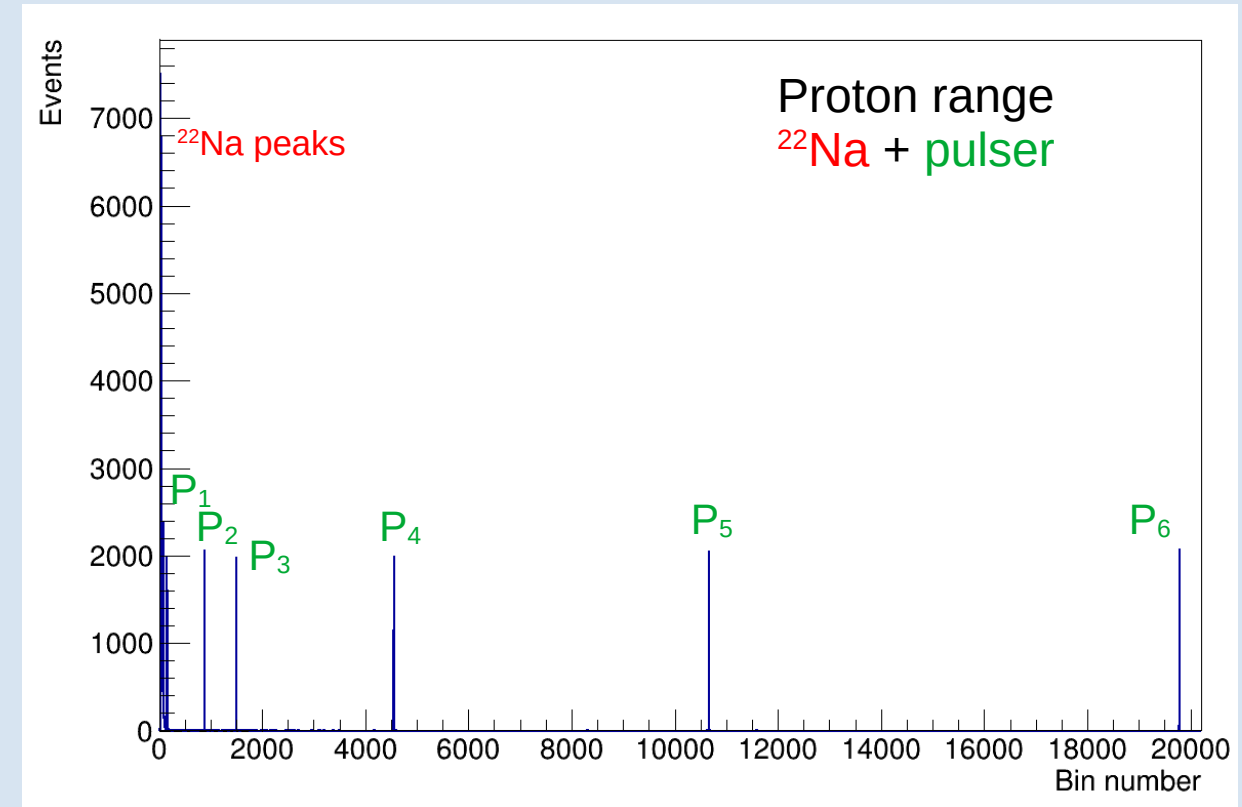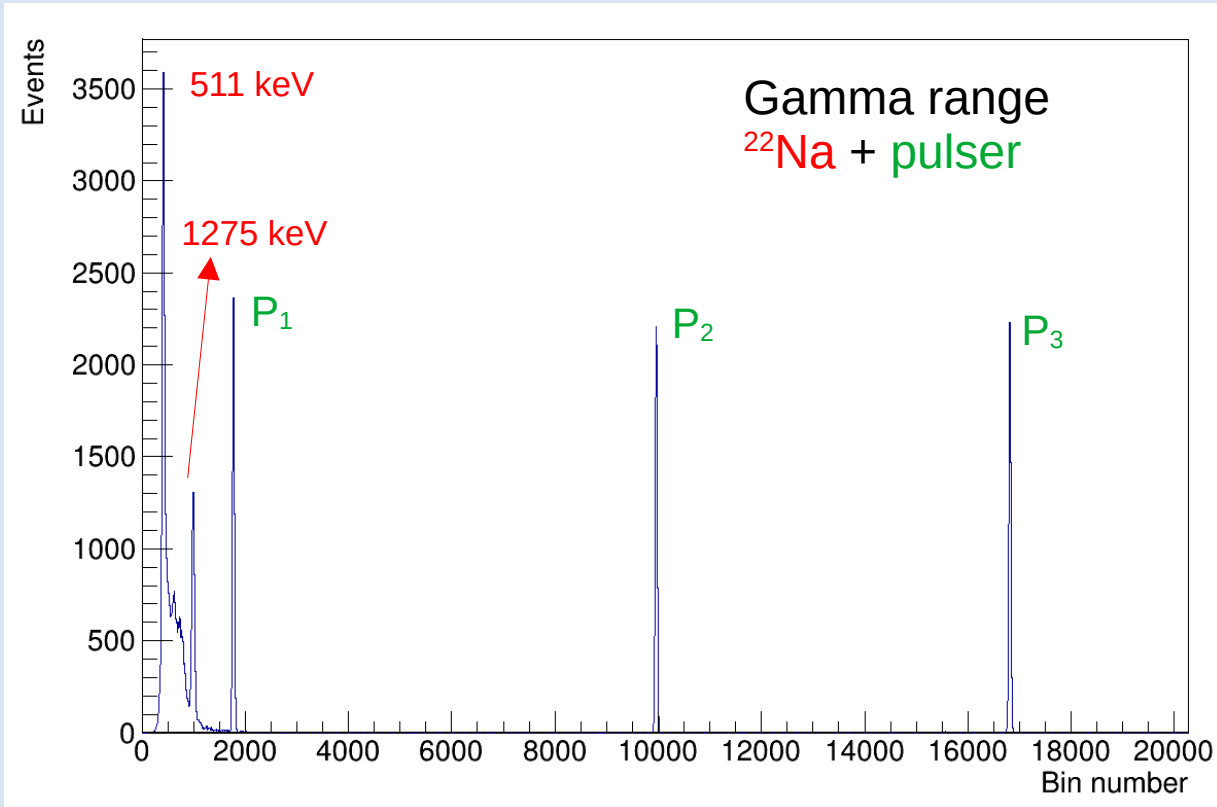

Using **pulsers** is essential for calibration in proton range

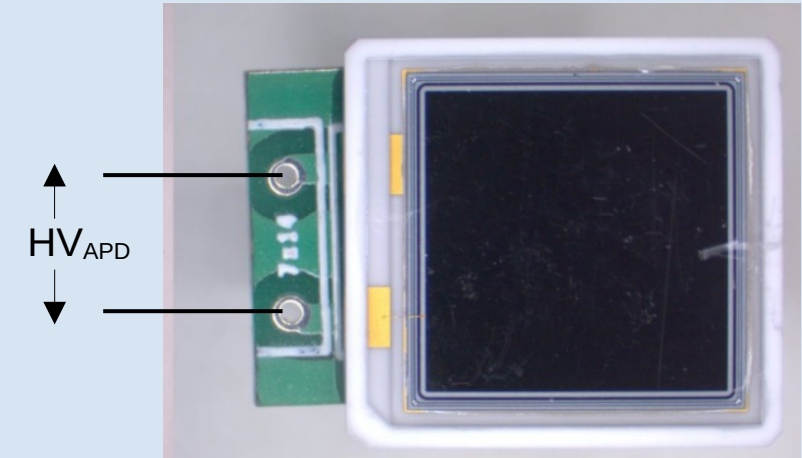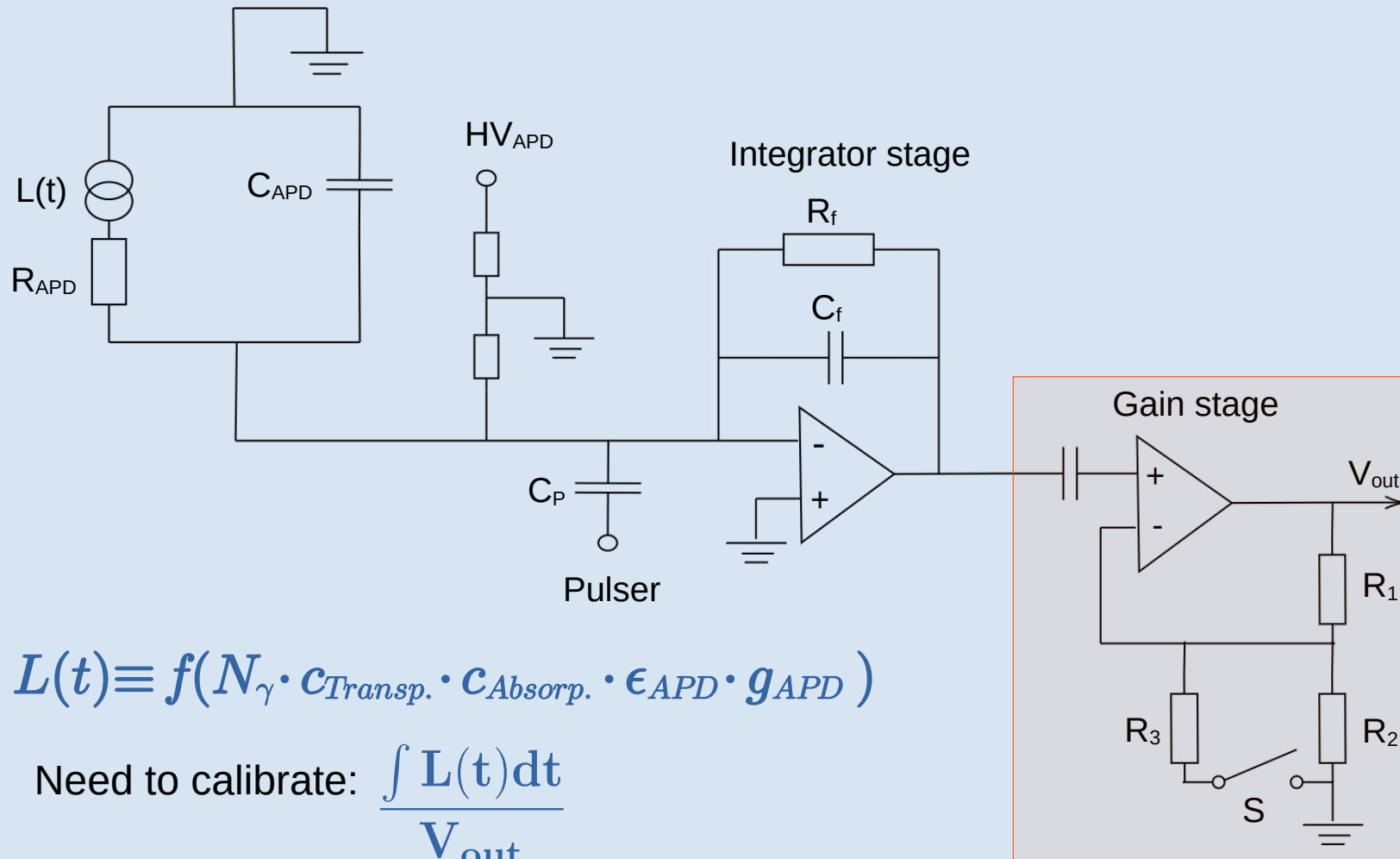


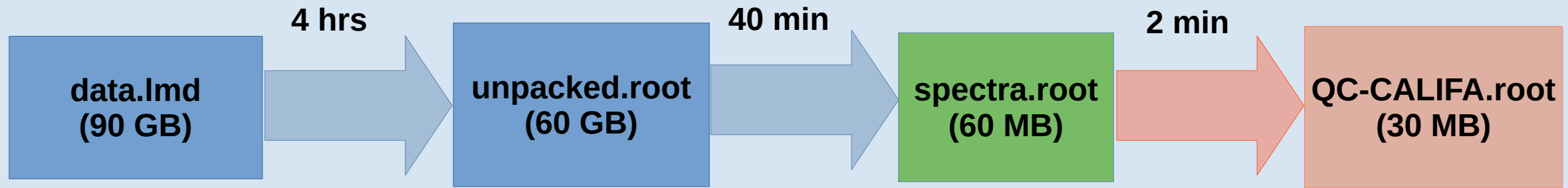

Key advantages:


- Unpacking, writing spectra done only once !
- Calibration takes just **2 min.**

Unpacking

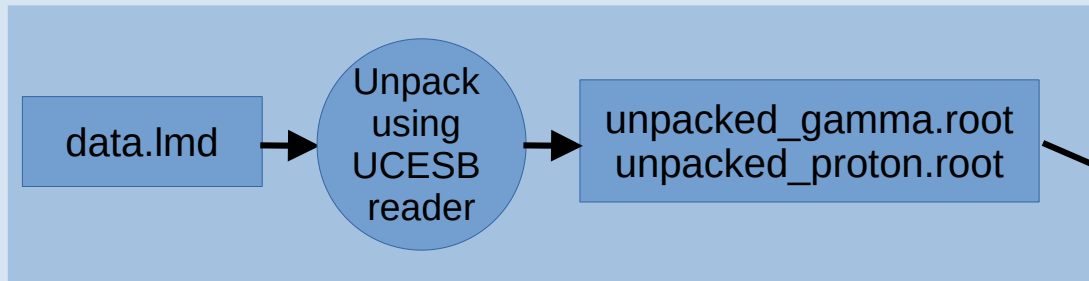


Source region


Pulser region

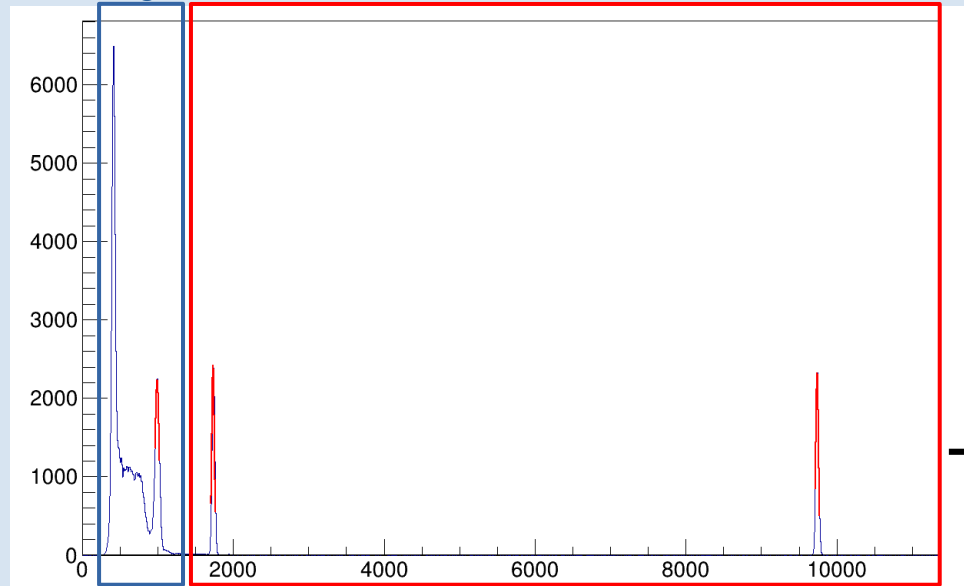



Calibration macro

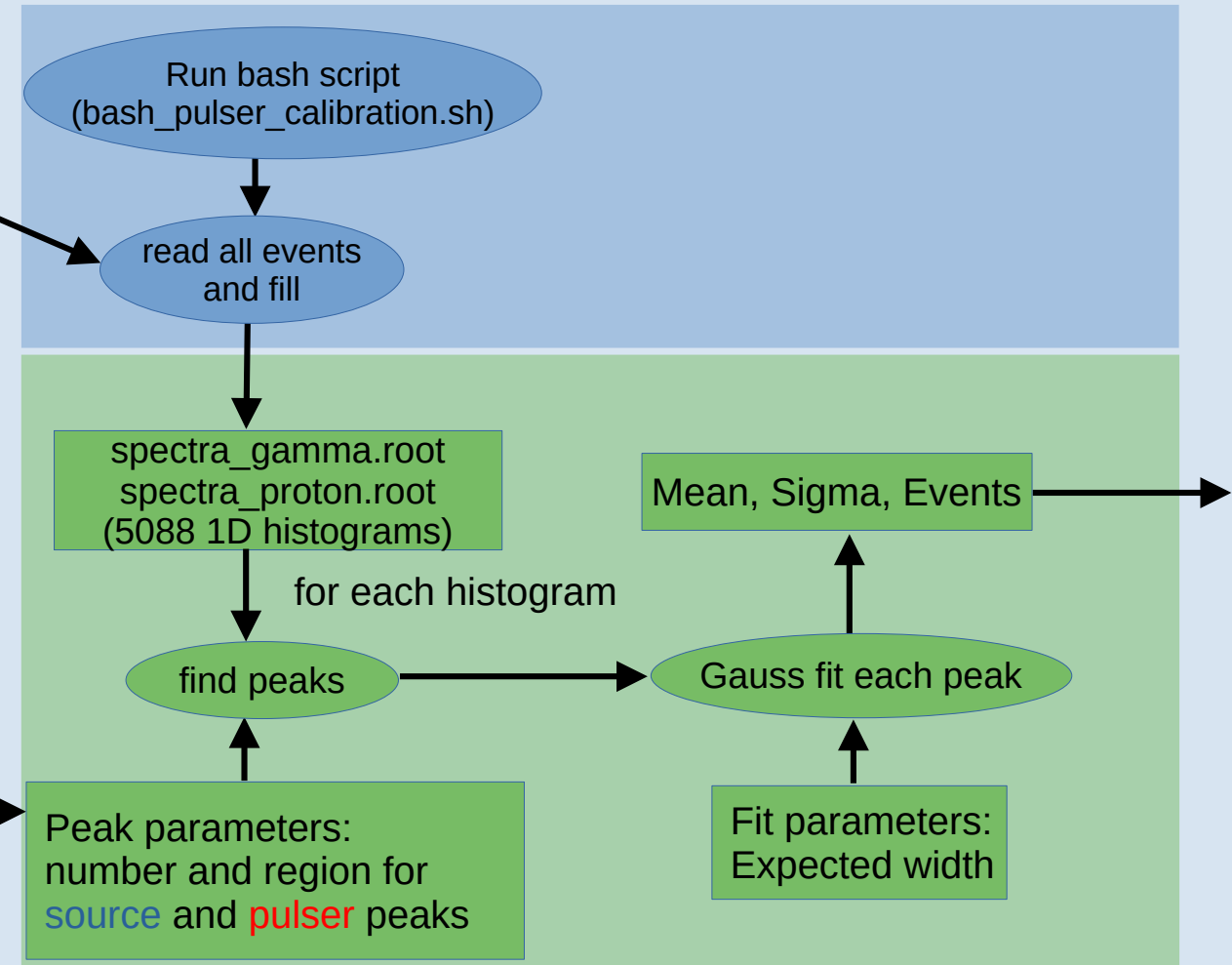

Calibration macro

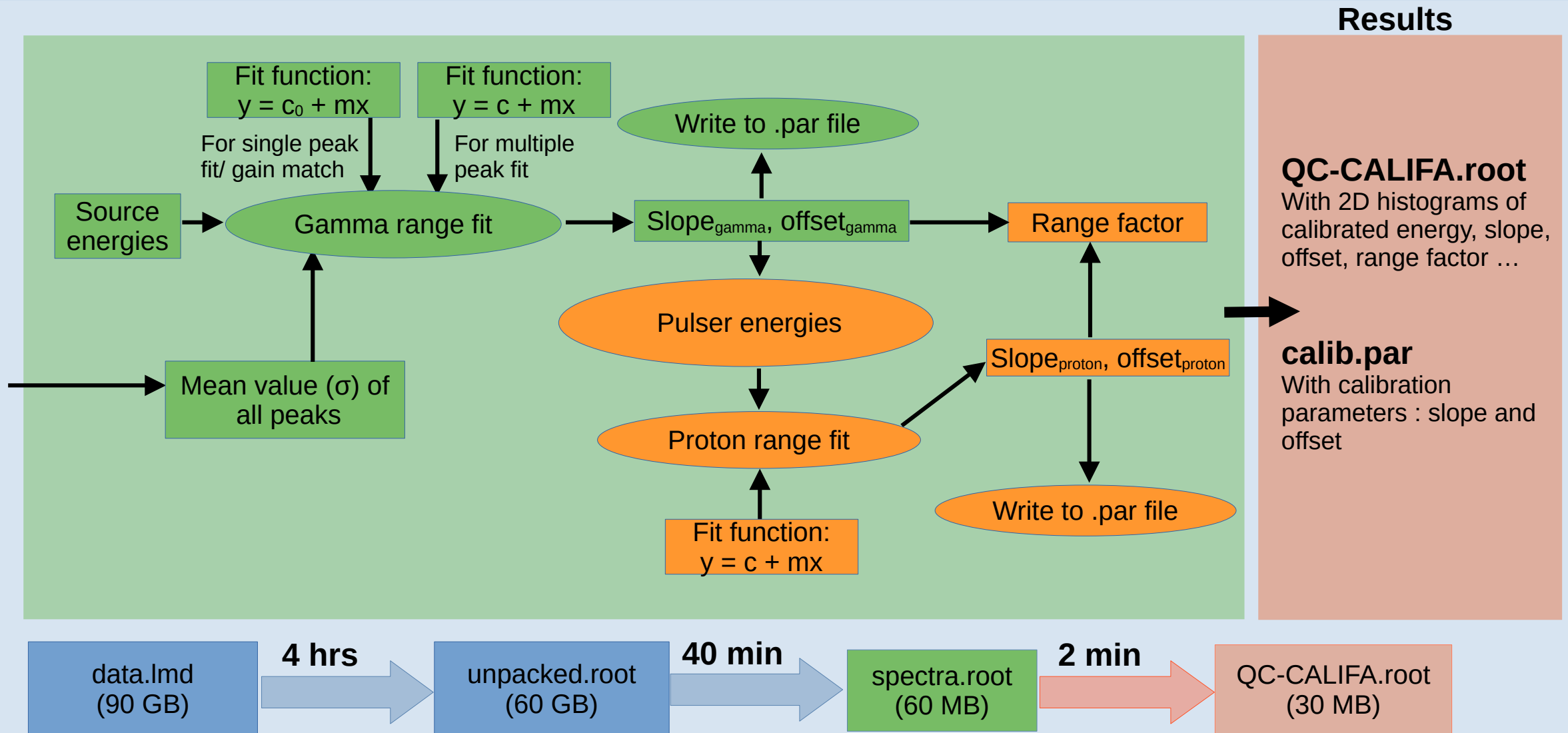




Range factor

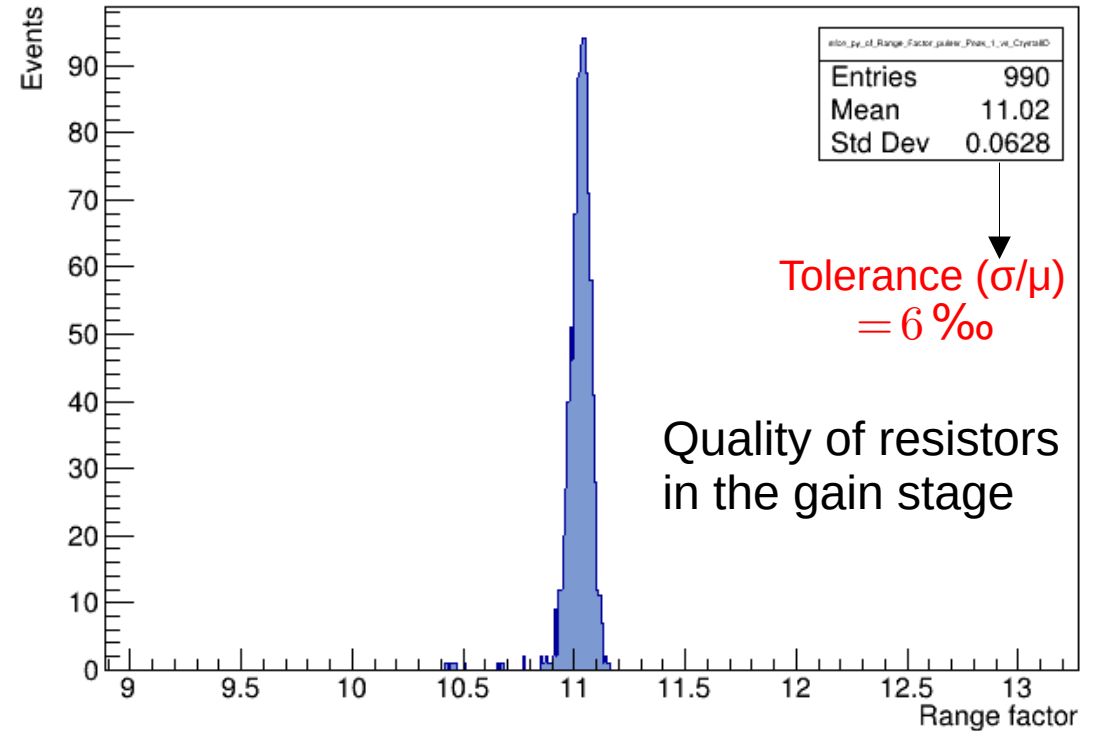
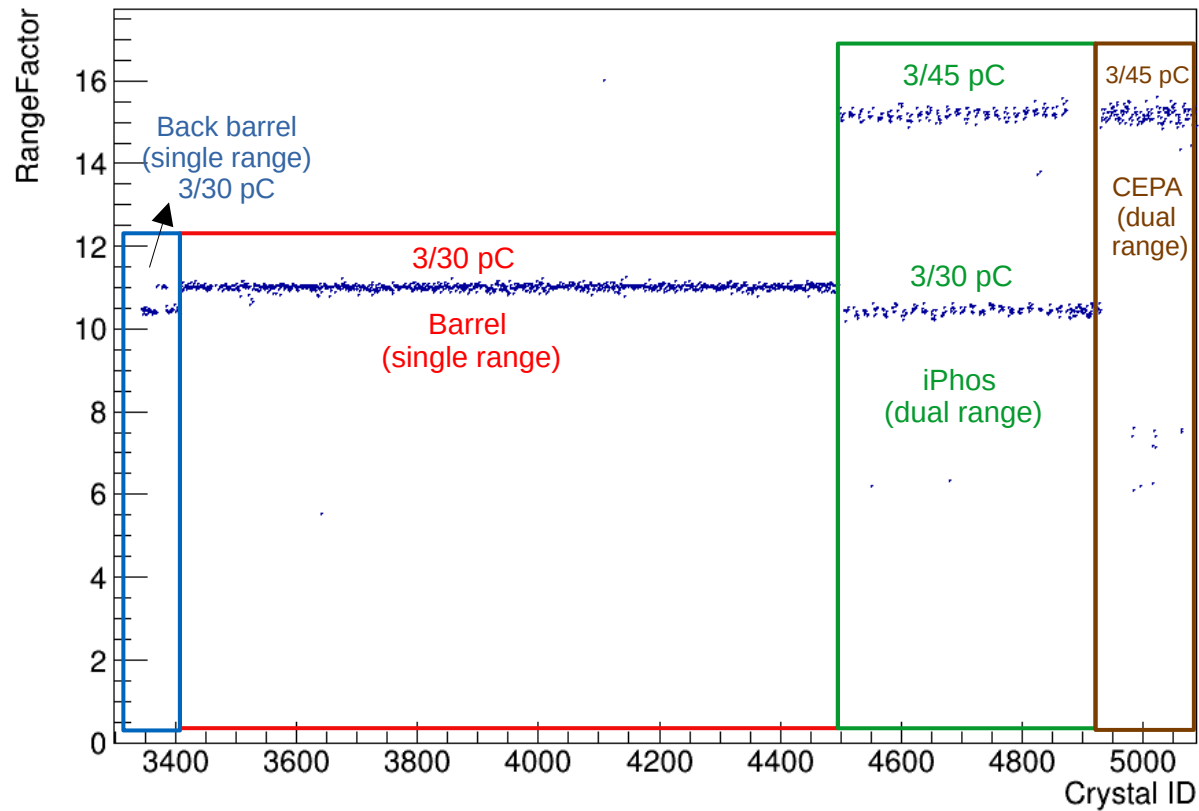


$$R = E_{\gamma(\text{uncalib.})} / E_{p(\text{uncalib.})} \approx 11$$

Serves as a translation factor between gamma range and proton range

Gamma range: uncalibrated spectra

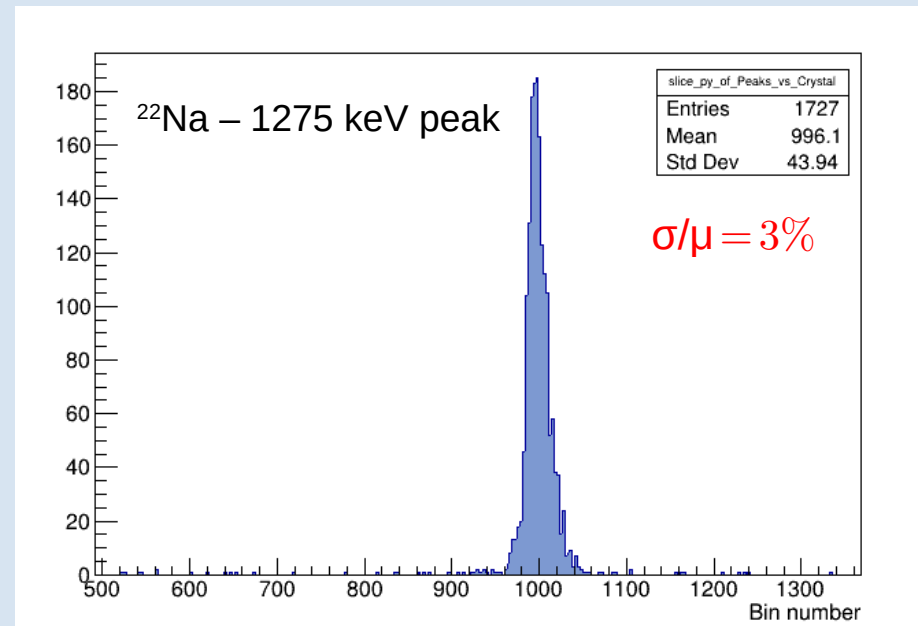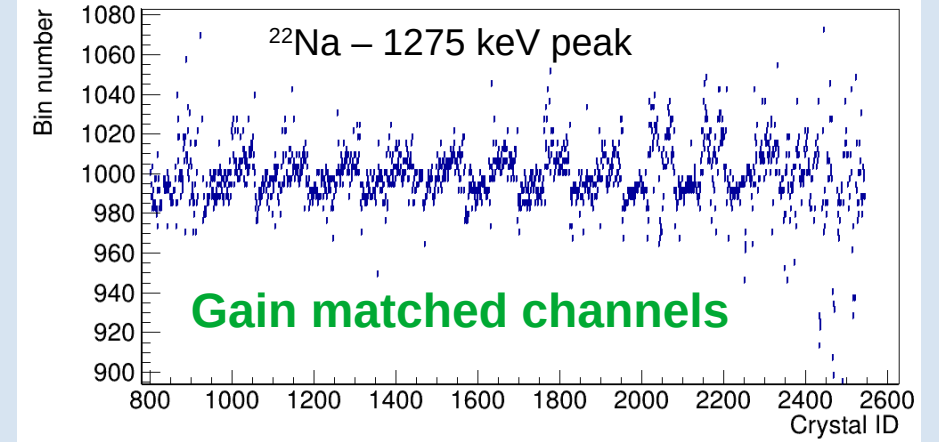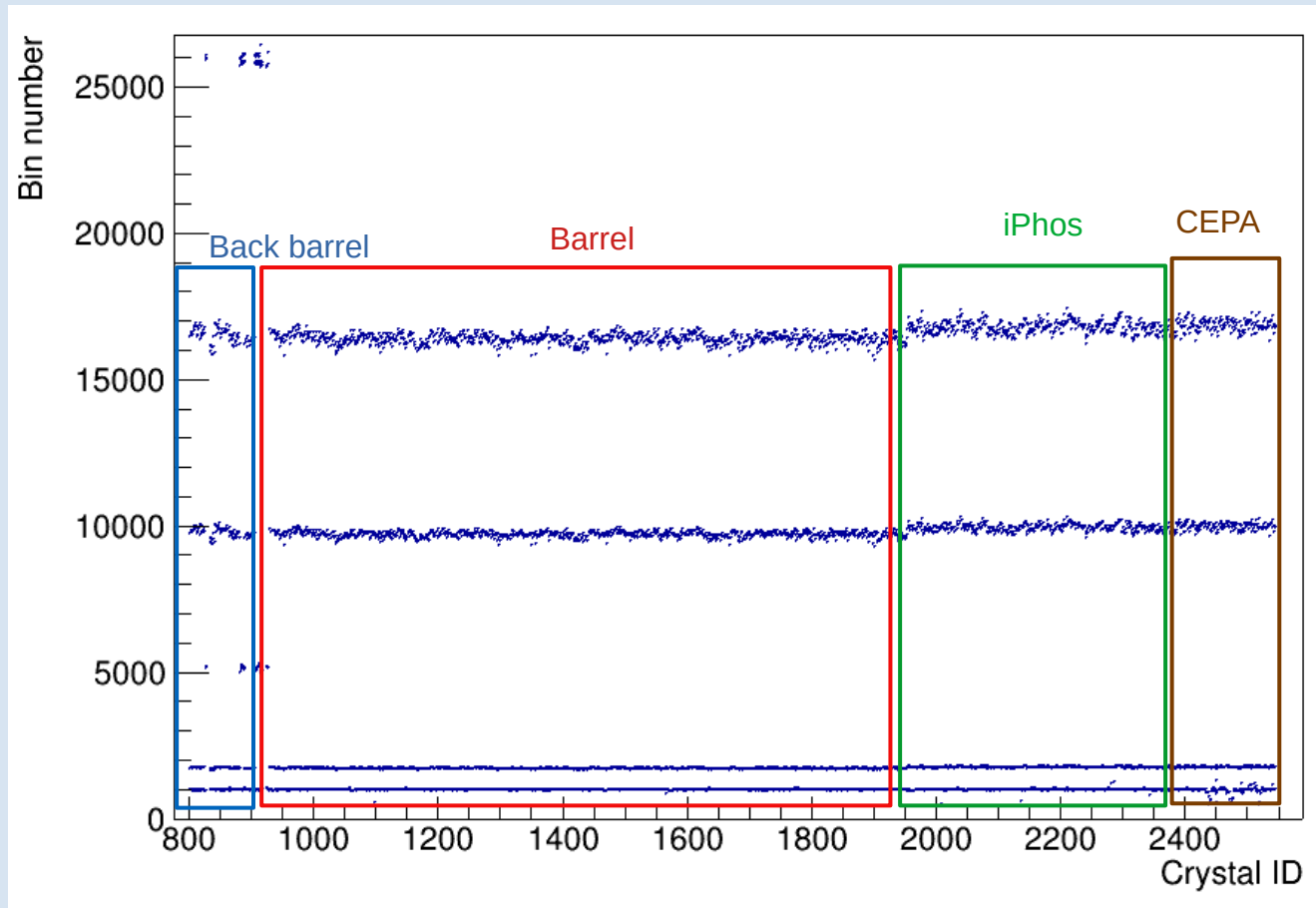

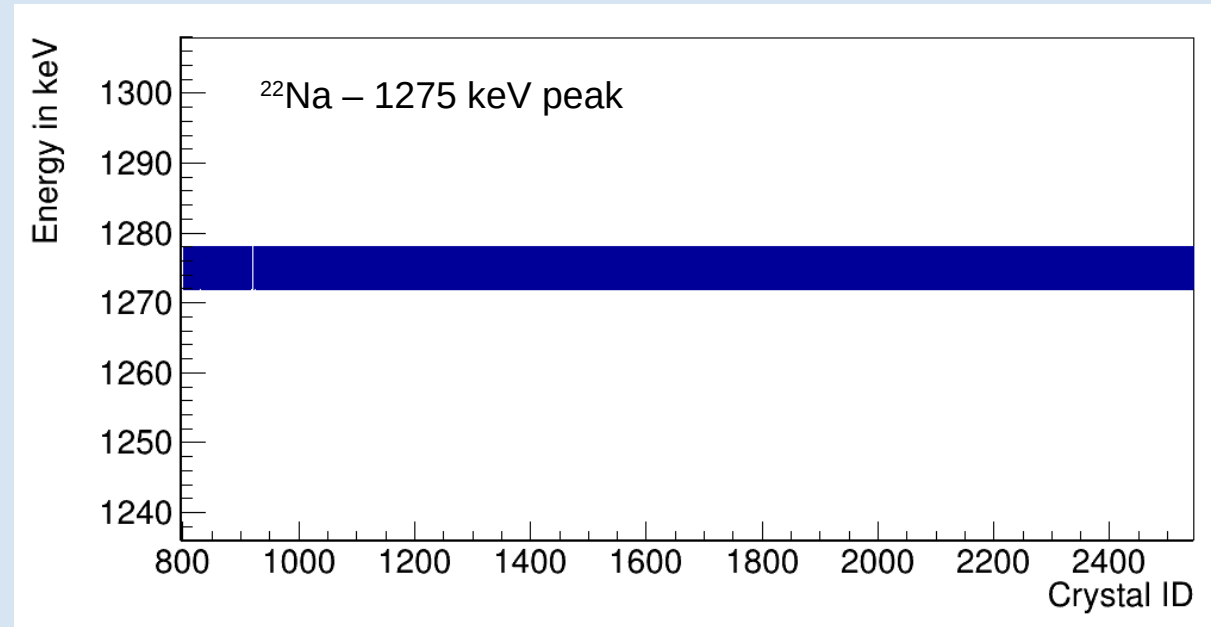
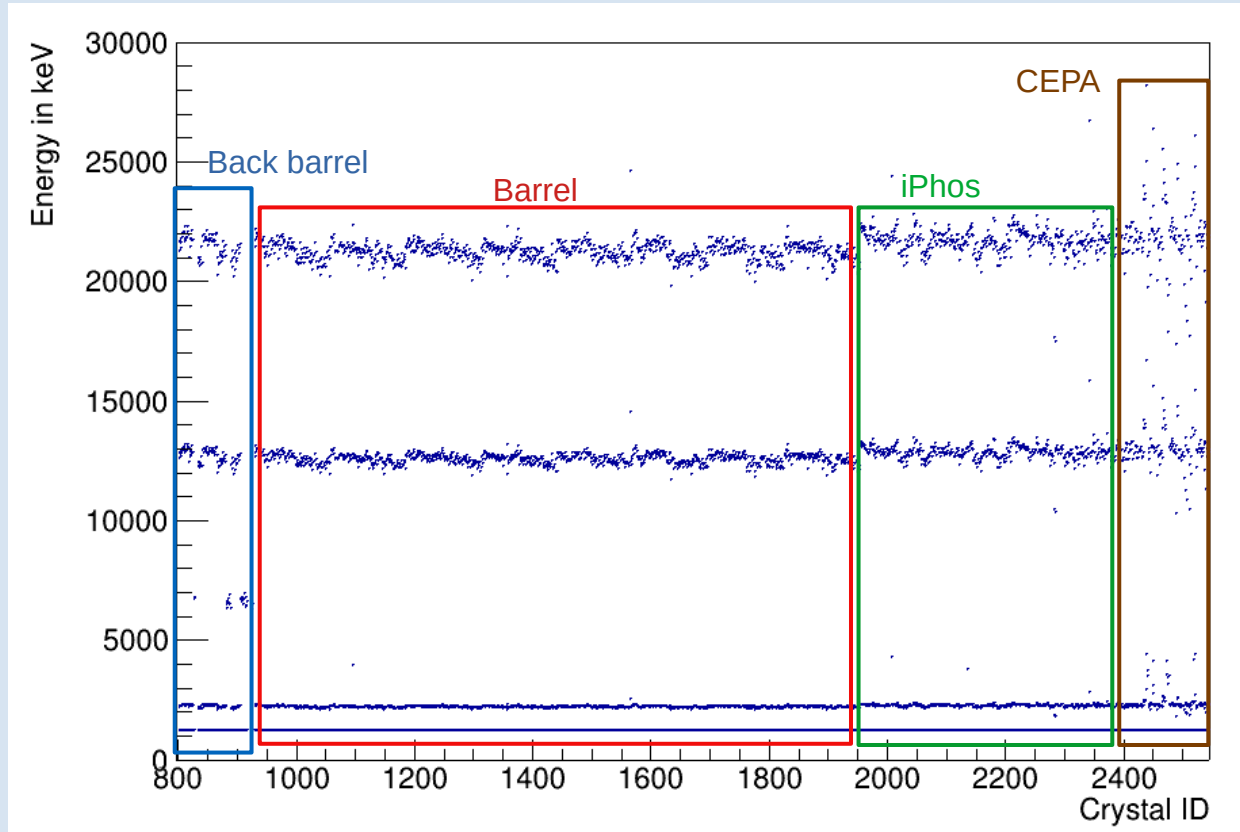


Applying a linear fit: $y = m_{\gamma}x + c_{\gamma}$



Software works !

Now we also have pulser peaks in keV

Proton range: uncalibrated spectra

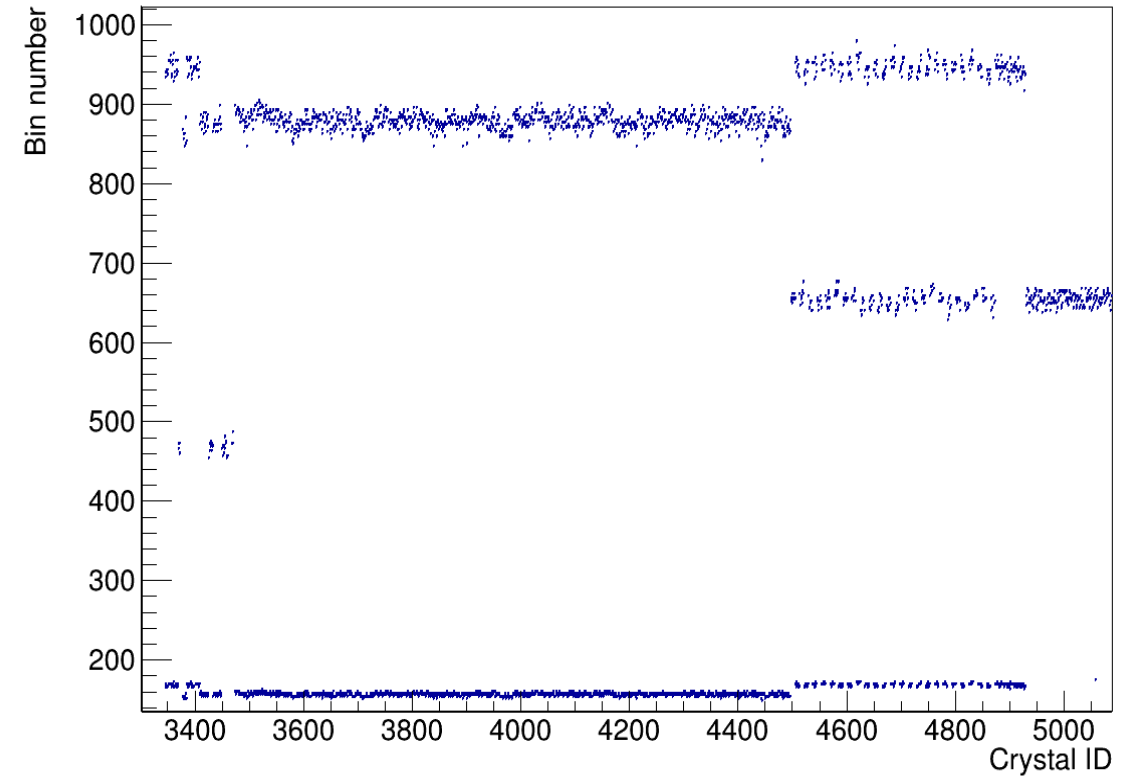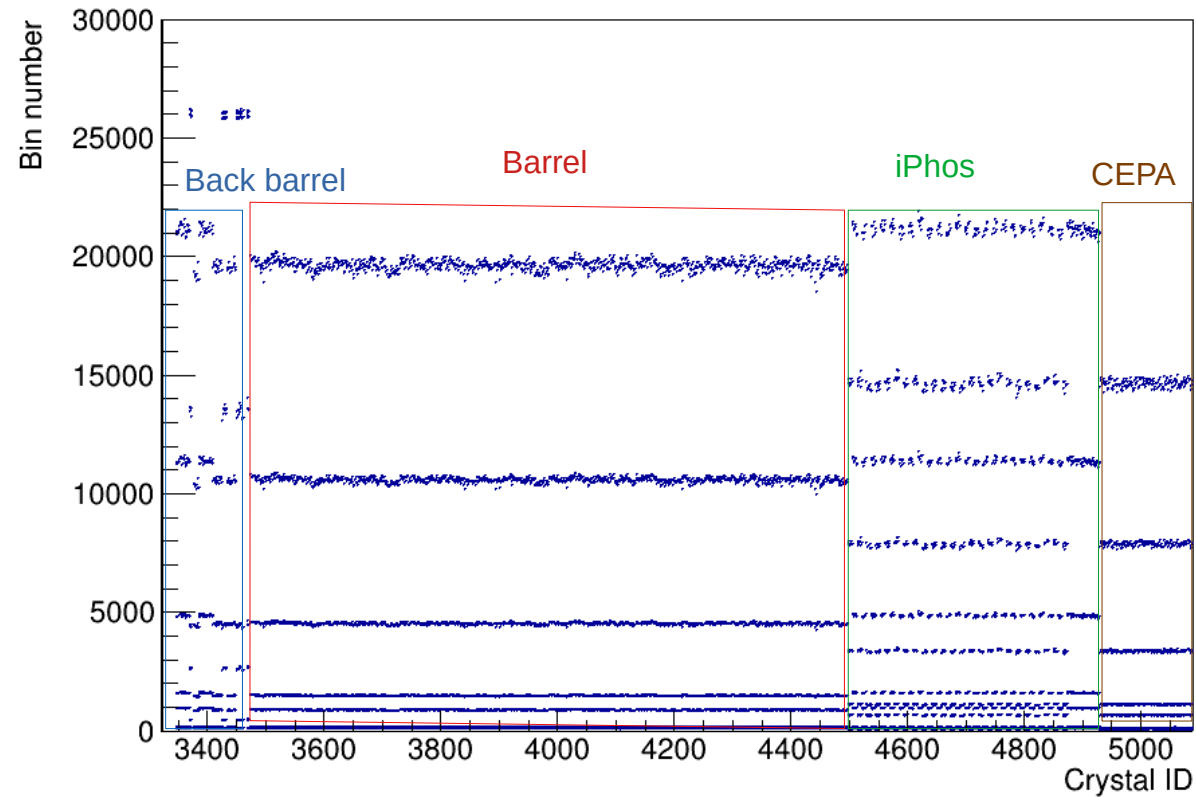


There are (channel to channel) variations in the electronics

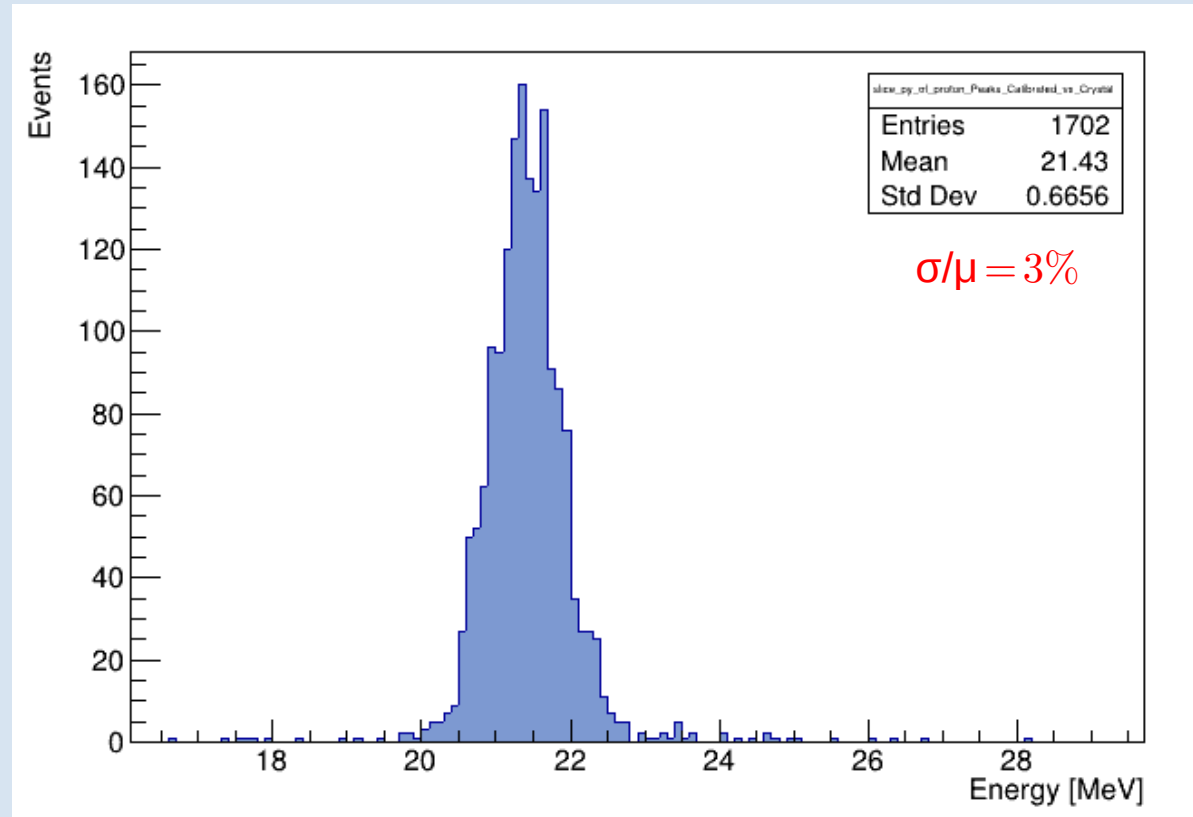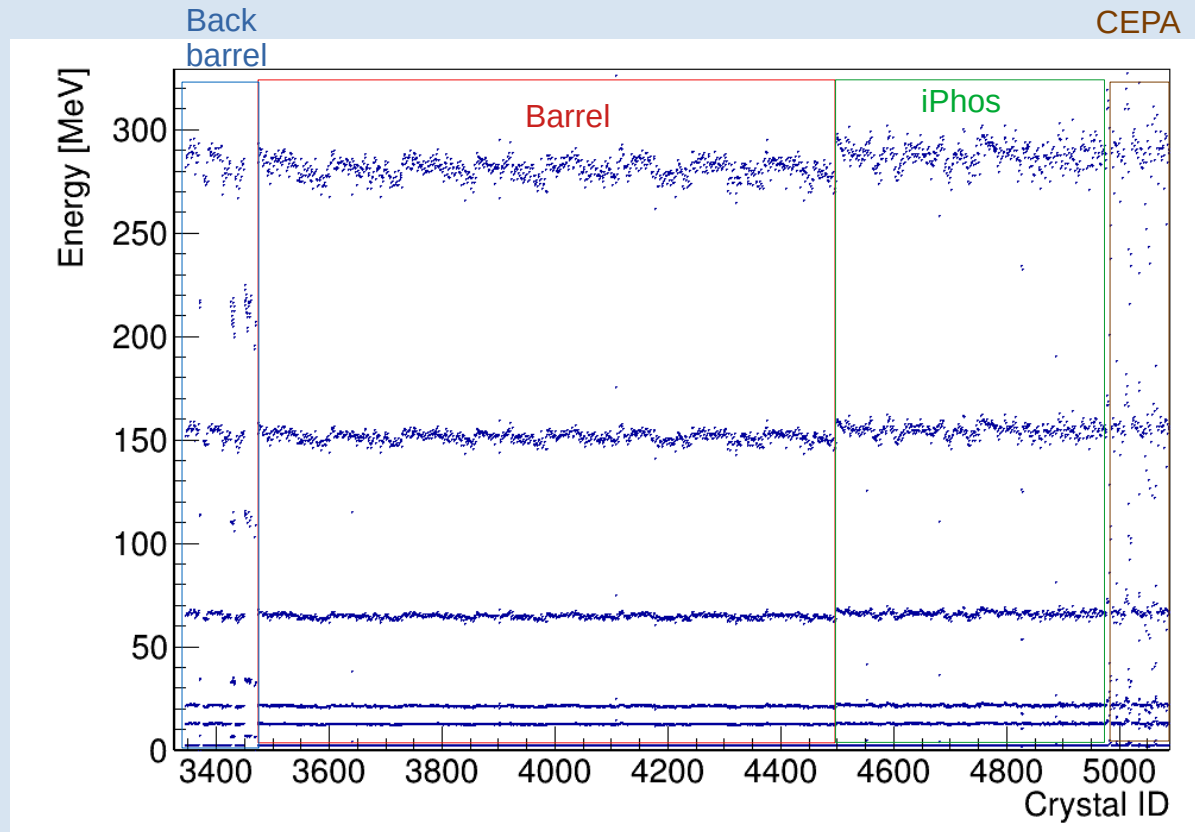


Energies for first 3 pulsers obtained from gamma calibration
Then doing a linear fit : $y = m_p x + c_p \Rightarrow$ All pulser energies in keV

Variation in pulsers include all uncertainties (depending on L(t), R and C_{pulser})

Circuit diagram


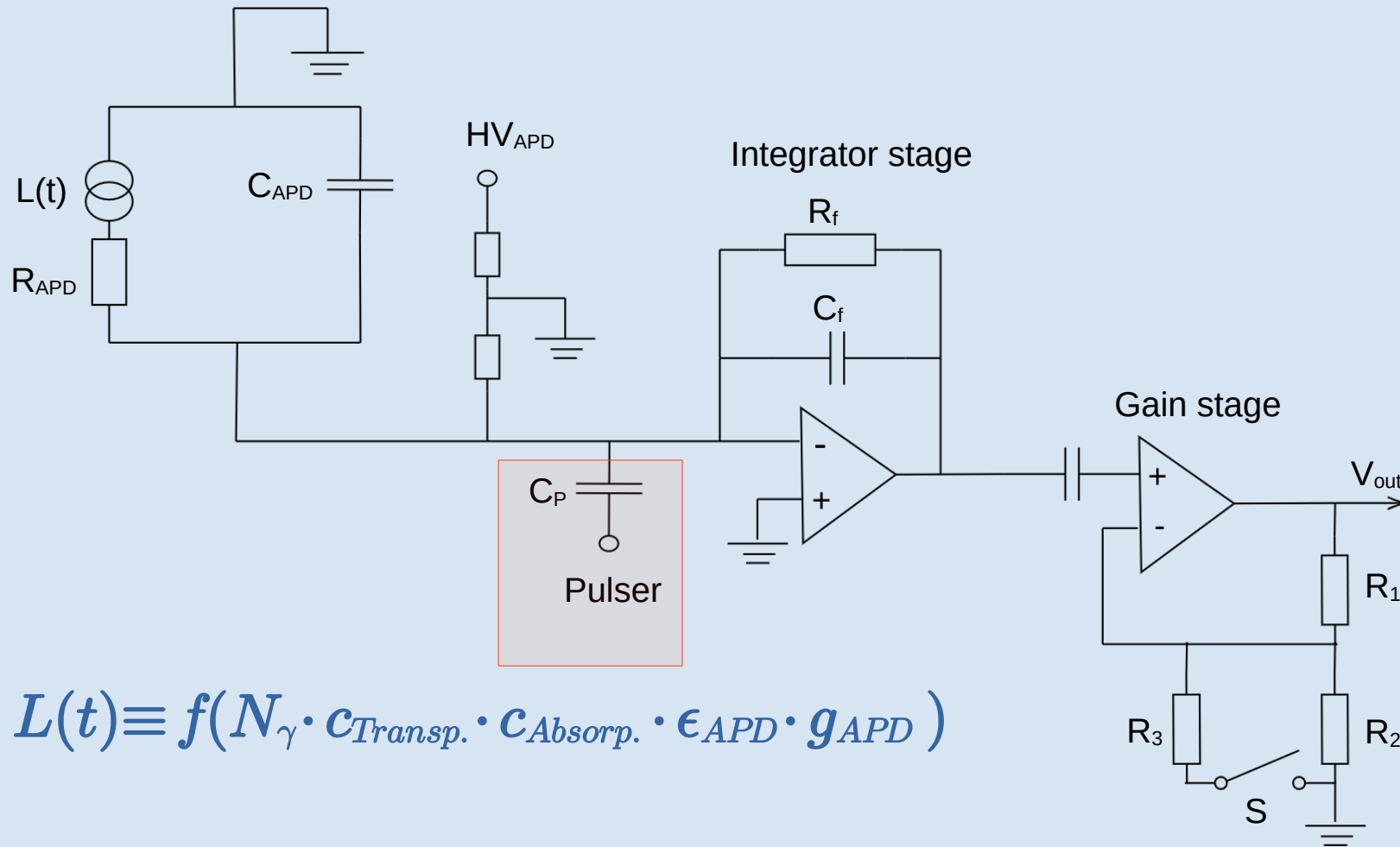


HV_{APD}

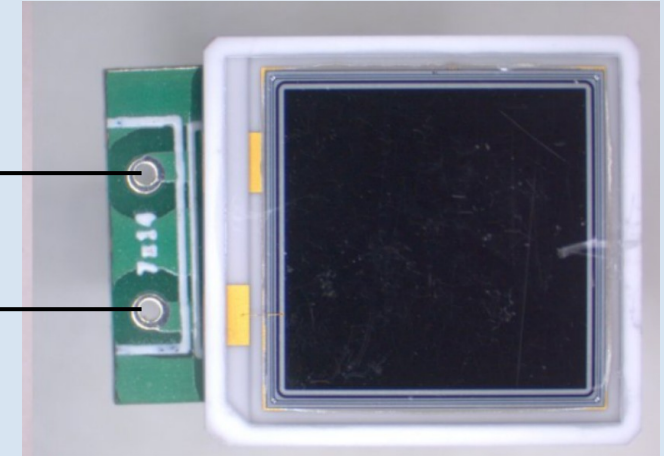



$$L(t) \equiv f(N_{\gamma} \cdot c_{Transp.} \cdot c_{Absorp.} \cdot \epsilon_{APD} \cdot g_{APD})$$

Pulser capacitance

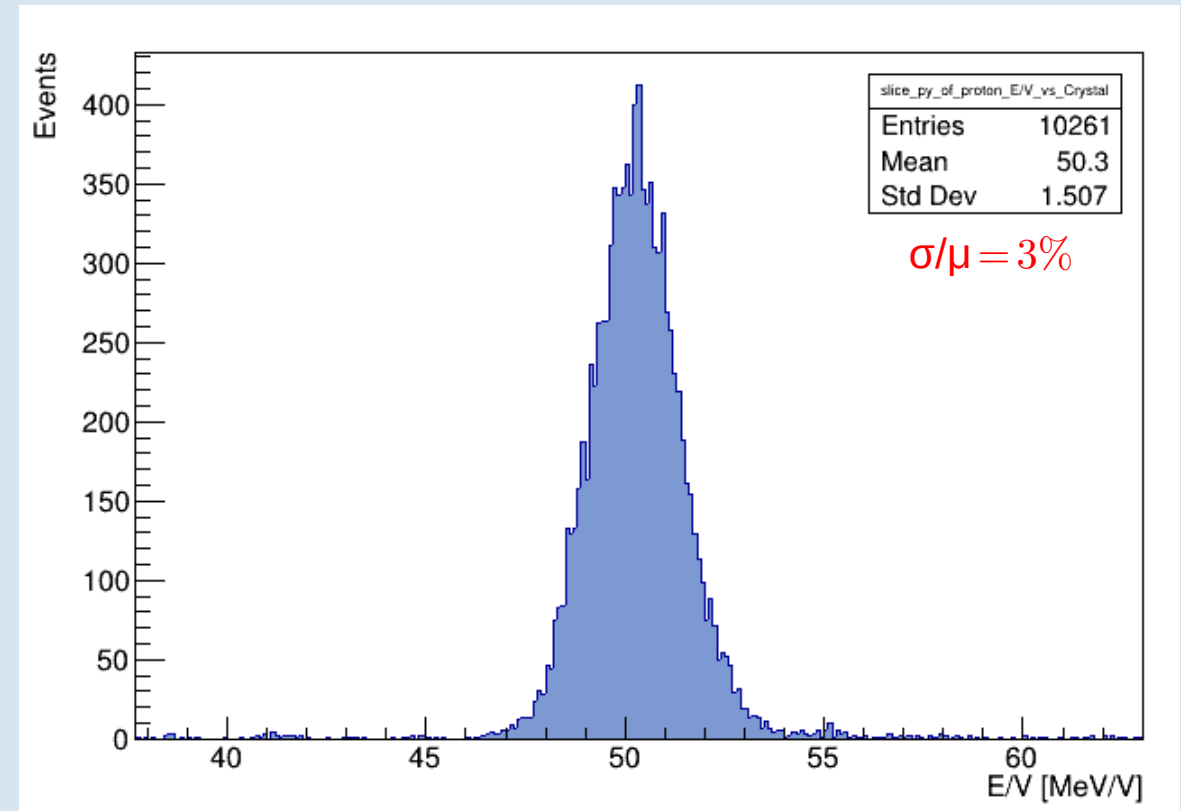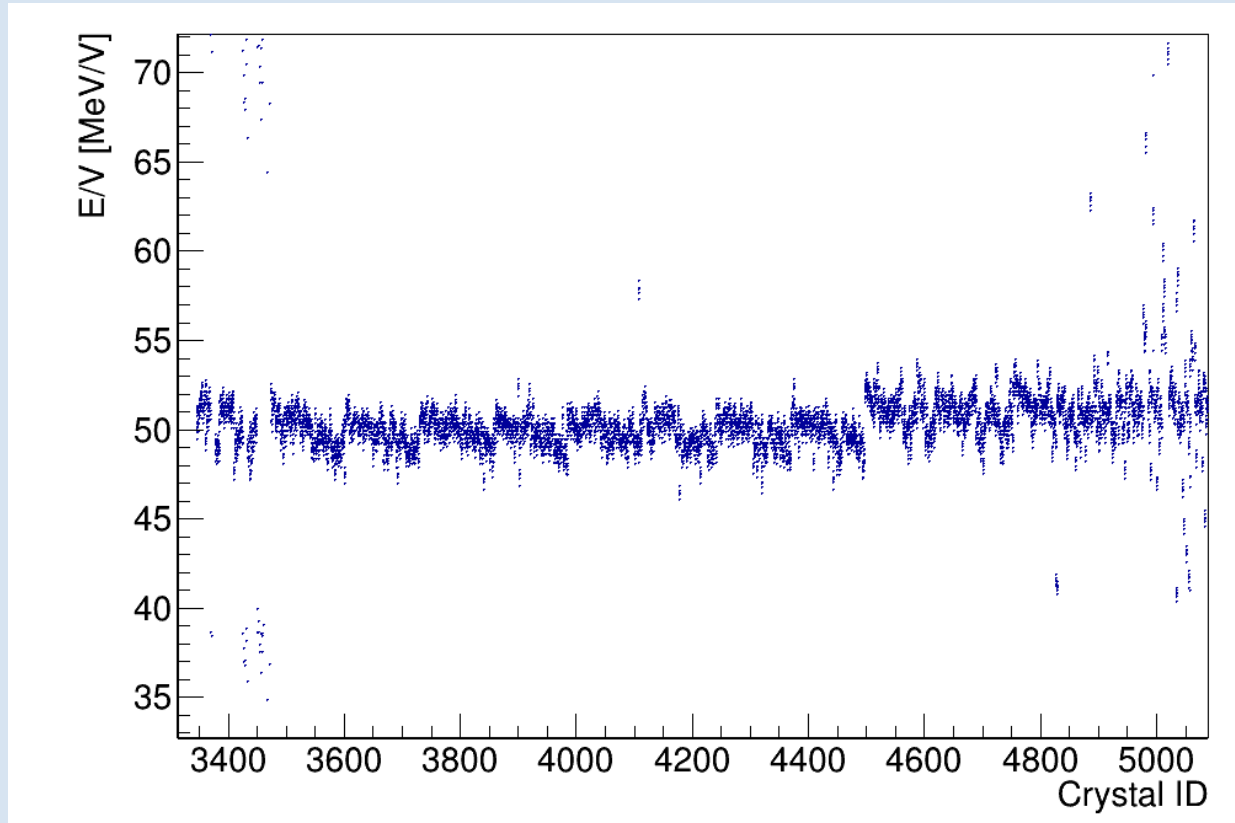


Variation in E/V for the pulsers again include all uncertainties (L(t), R and C_{pulser})


Proton range: calibrated spectra

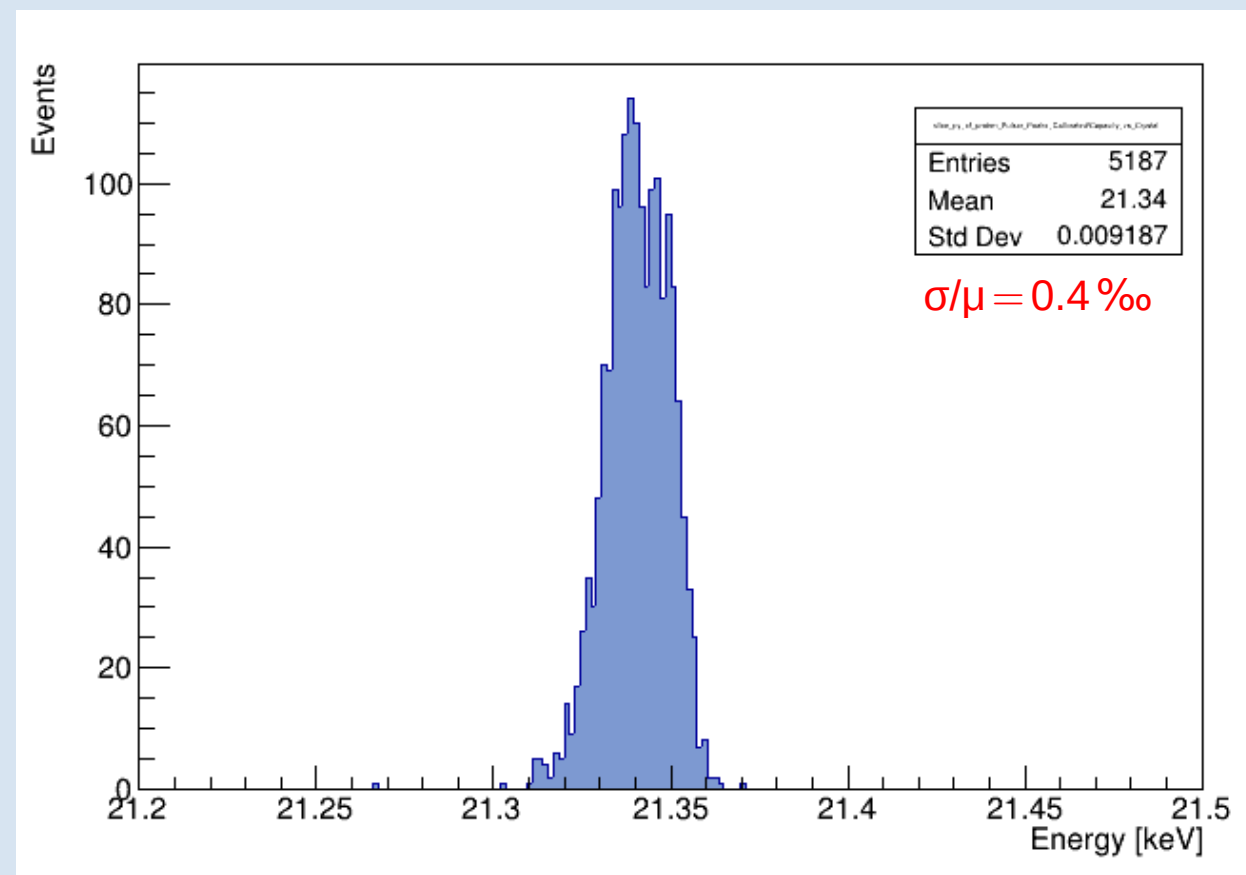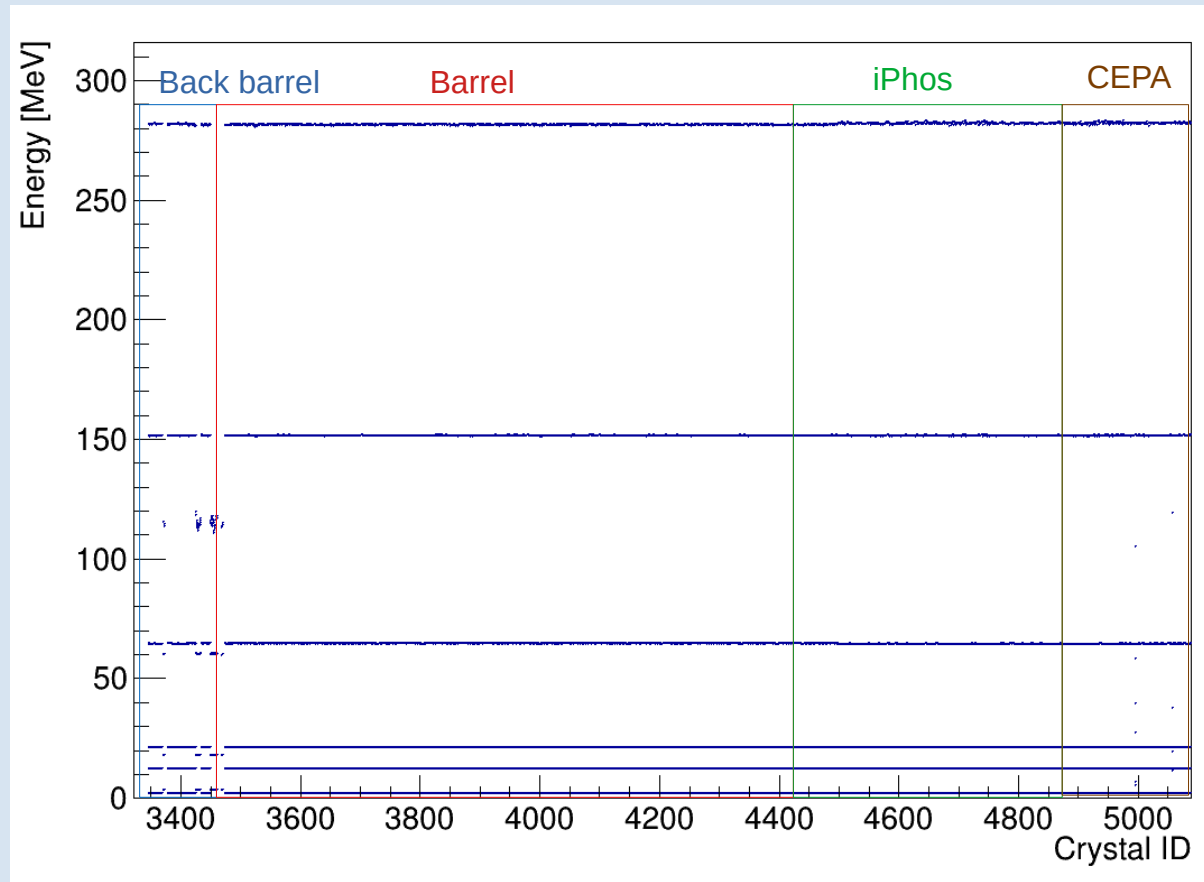



Software works !

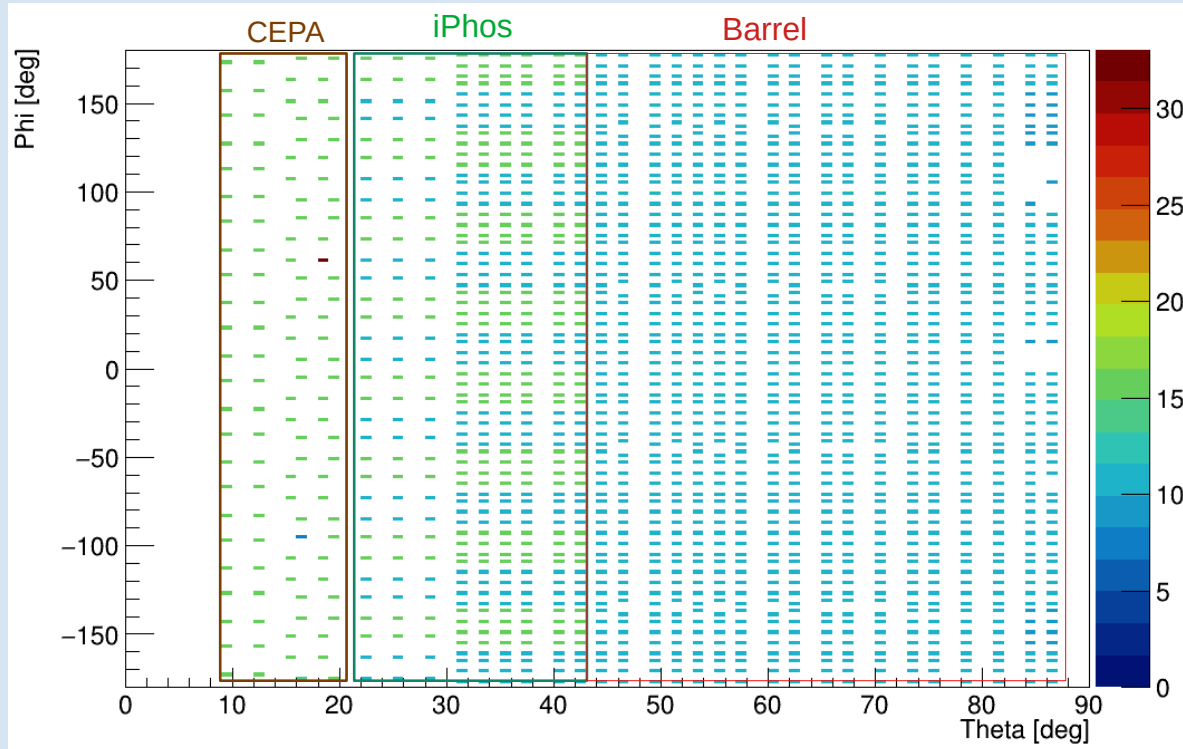

Quality Assessment Plot

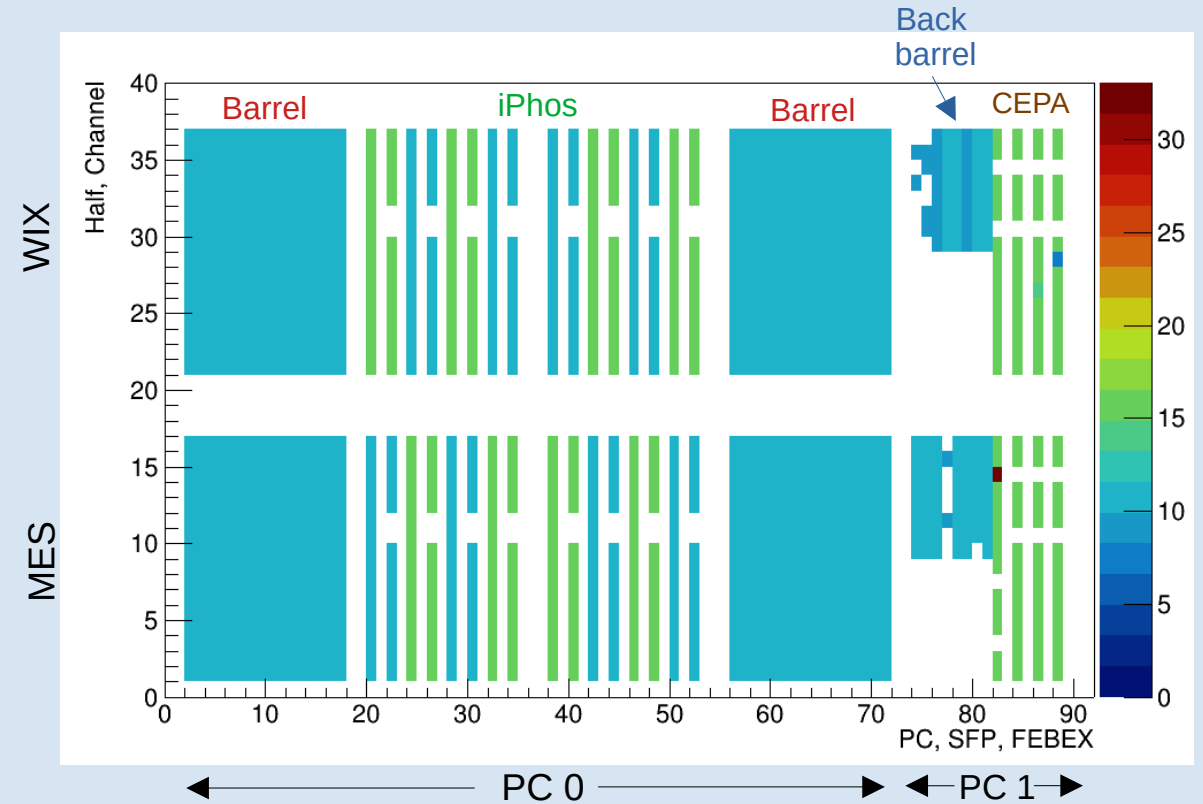

DAQ Plot (Hardware level)

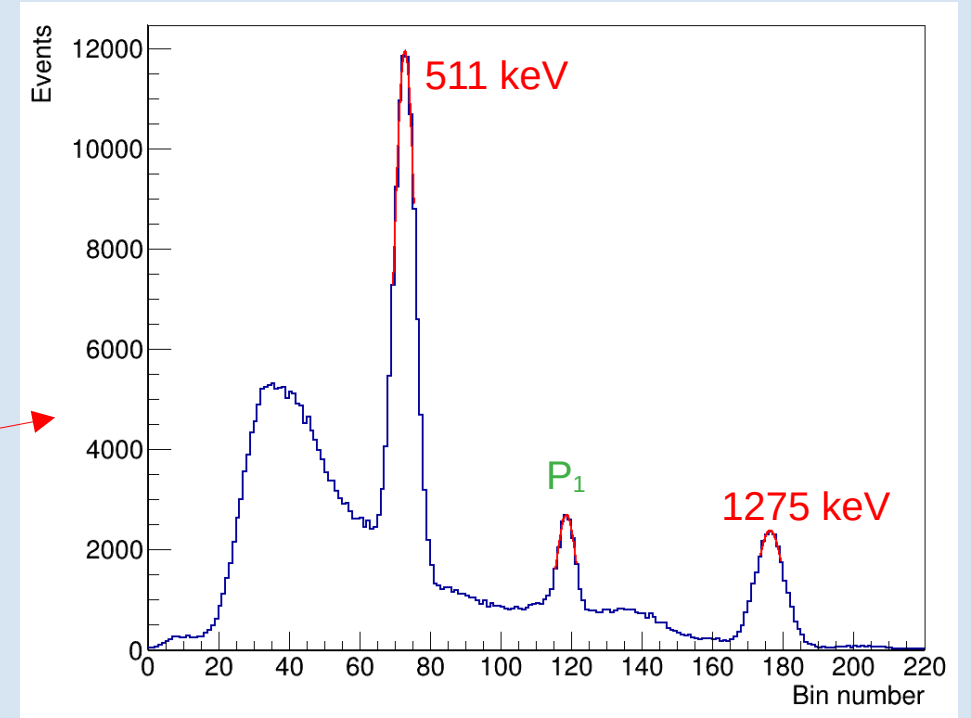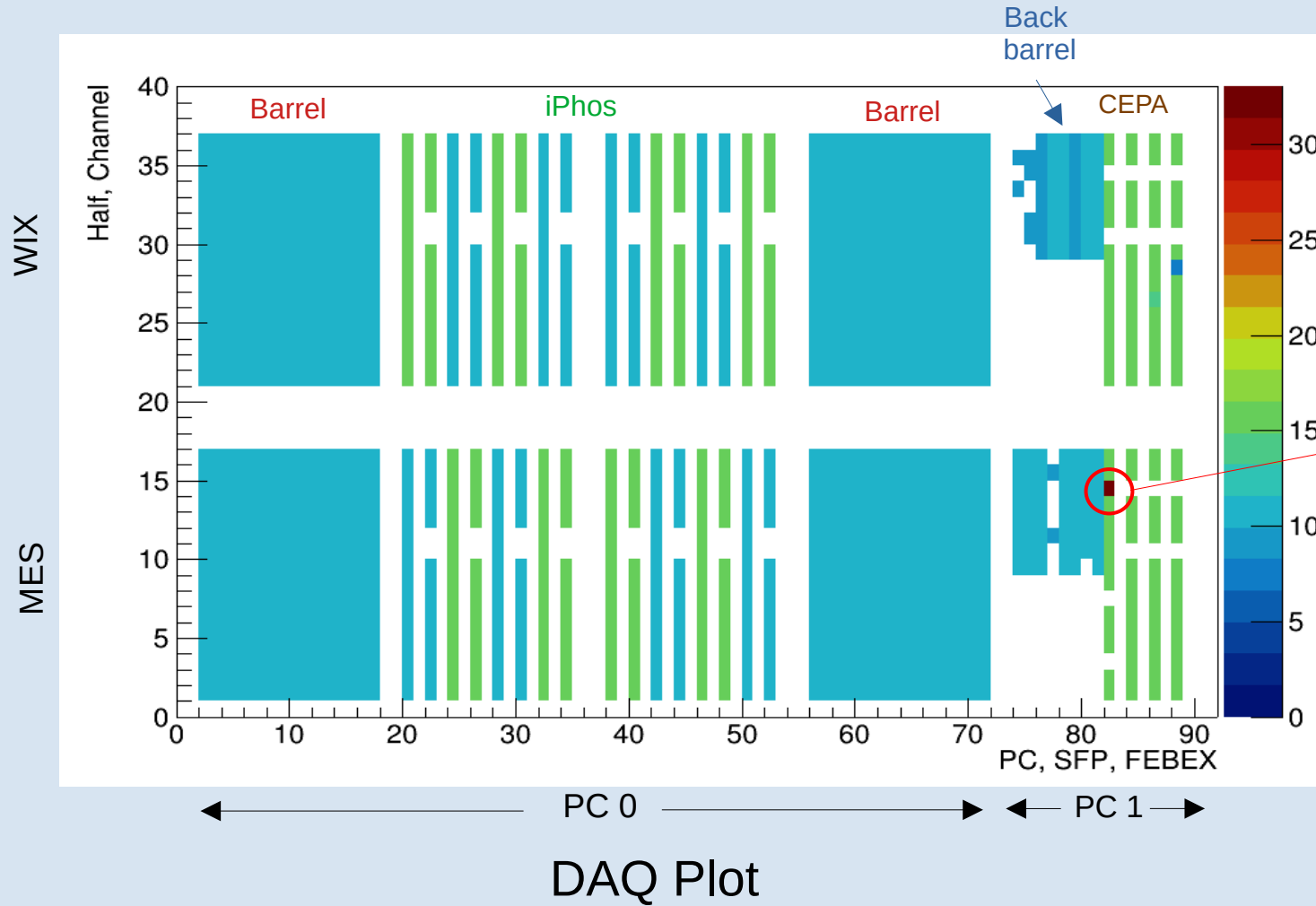

Spectrum for an incorrectly gain matched channel




Thank You !

CALIFA @ Technical University of Munich

Roman Gernhäuser, Philipp Klenze, Tobias Jenegger, Mrunmoy Jena



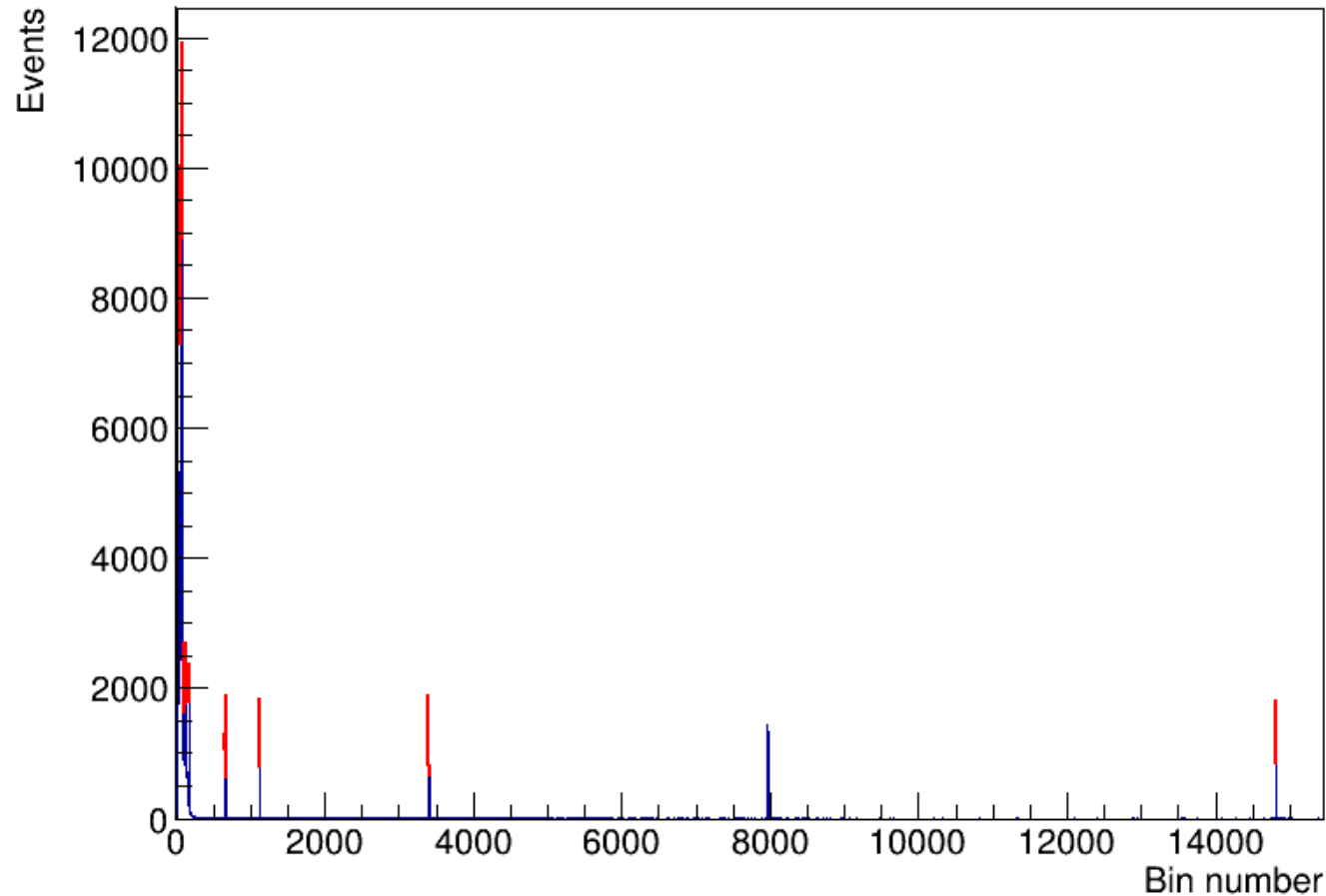

Crystal ID: 5057, in CEPA

Crystal 5057

Source Peak 1, Bin number: 72.8852
Pulser Peak 1, Bin number: 118.622
Pulser Peak 2, Bin number: 176.335
Pulser Peak 3, Bin number: 664.026
Pulser Peak 4, Bin number: 1121.33
Pulser Peak 5, Bin number: 3402.16
Pulser Peak 6, Bin number: 14819.5
Range factor Pulser Peak 1: 15.2875
Range factor Pulser Peak 2: 57.8574
Range factor Pulser Peak 3: 25.9151
Range factor: 33.02
Pulser Offset: -2053.82
Pulser Slope: 36.794

Extras: Range factor vs bin number

